

3M Science.
Applied to Life.™

3M™ ESPE™ Sof-Lex™ Diamond Polishing System

Technical Data Sheet





Product Description

The **3M™ ESPE™ Sof-Lex™ Diamond Polishing System** is a two-step, multi-use polishing system. The spirals use the same mandrel as 3M™ ESPE™ Sof-Lex™ Finishing and Polishing Discs. The spirals are made with either aluminum oxide or diamond particles impregnated in a thermoplastic elastomer.

The universal shape allows for usage on all tooth surfaces and reduces the need for different shaped tools (eg. points, cups, discs and brushes), which are designed for specific shapes and contours.

One shape:

- Adapts to all tooth surfaces
- Works from any angle
- Effective for anterior and posterior restorations
- Helps create natural gloss in two steps

Indications

The **3M™ ESPE™ Sof-Lex™ Pre-Polishing Spiral** (beige) is indicated to smooth the surface of restorations. The **3M™ ESPE™ Sof-Lex™ Diamond Polishing Spiral** (pink) is indicated for final high-gloss polishing. Spirals can be used to polish surfaces of:

Direct:

- Composite restorations
- Resin-modified glass ionomers
- Bisacrylic temporary materials

Indirect:

- Composite
- 3M™ ESPE™ Lava™ Ultimate CAD/CAM Restorative
- Precious and semi-precious metal

Both wheels are used with a slow-speed handpiece operating within 15,000–20,000 rpm. Adding water is recommended during use and will create a more ideal surface quicker than without water.


Introduction:

Proper finishing of restorations is desirable not only for esthetic considerations but also for oral health. The primary goal of finishing is to obtain a restoration which has ideal contour, occlusion, healthy embrasure forms and surface smoothness. Finishing and polishing procedures remove the air-inhibited layer, contour and shape the restoration, create surface characterization and produce surfaces with high gloss. Tight margins should blend esthetically into the tooth's natural contours. Healthy embrasure forms and smooth surfaces are less likely to trap food debris and collect plaque. Several factors can affect the final finish of a restoration: the resin matrix and fillers within the restorative material, finishing instruments and preparation design.

General Procedure

A full line of finishing and polishing solutions ... from beginning to end.


The Sof-Lex system is color coded from dark (coarse) to light (superfine) for an easily followed step-by-step process. Easy, pop-on mounting lets you quickly change discs or spirals—so you can move efficiently through your finishing and polishing procedure.

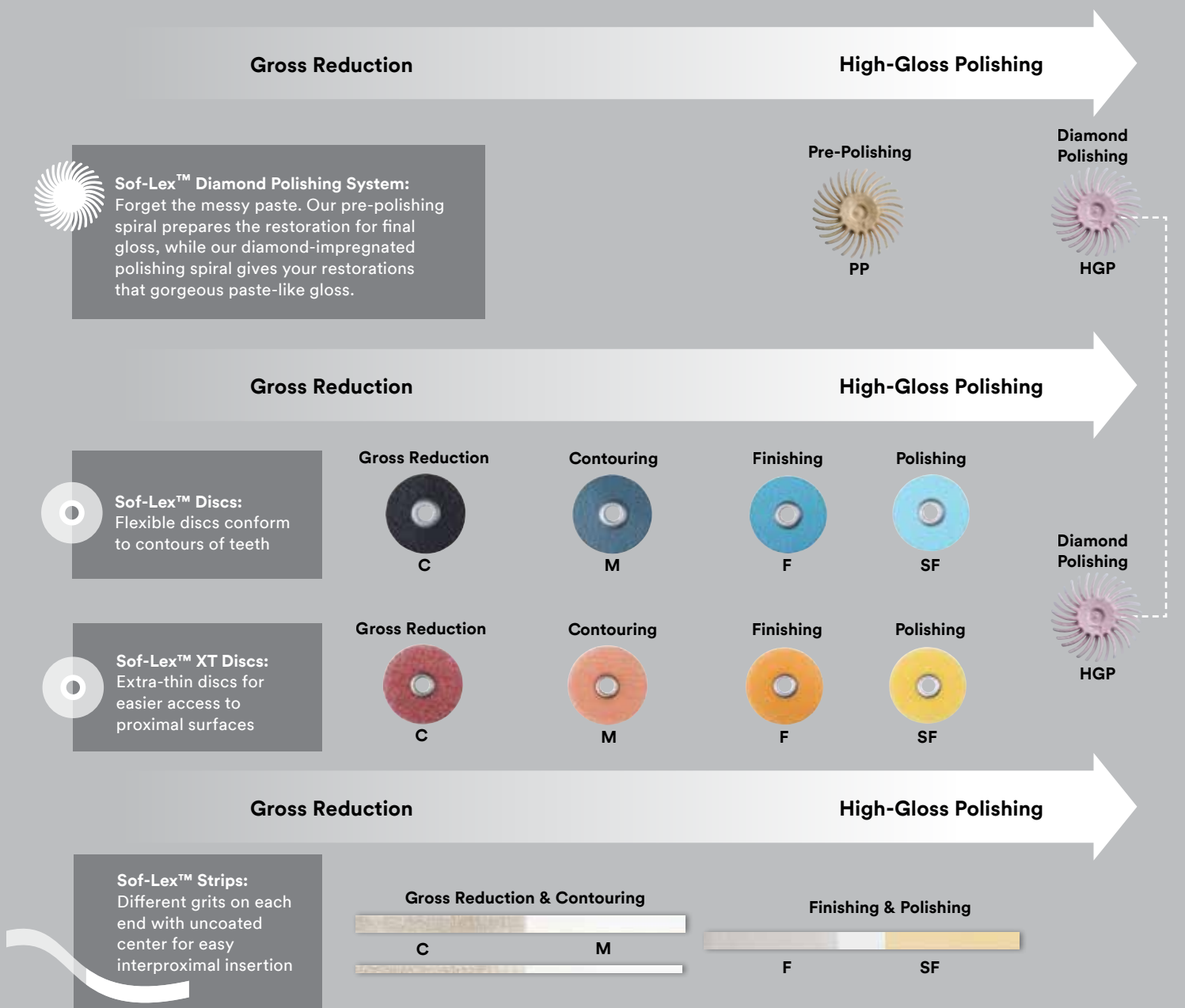

Gross Reduction

The purpose of the Gross Reduction step is to remove excess restorative material (including overhangs), to remove the air-inhibited layer and shape anatomy. The tools that are typically used for this process deliver the most aggressive abrasive (e.g. coarse grit) action in this procedure to remove the excess material quickly. The graph below depicts the most popular tools for Gross Reduction. Abrasive strips are also used for proximal areas.

Percentage of respondents who use the corresponding tools for the Gross Reduction step (N=300). Multiple responses allowed.

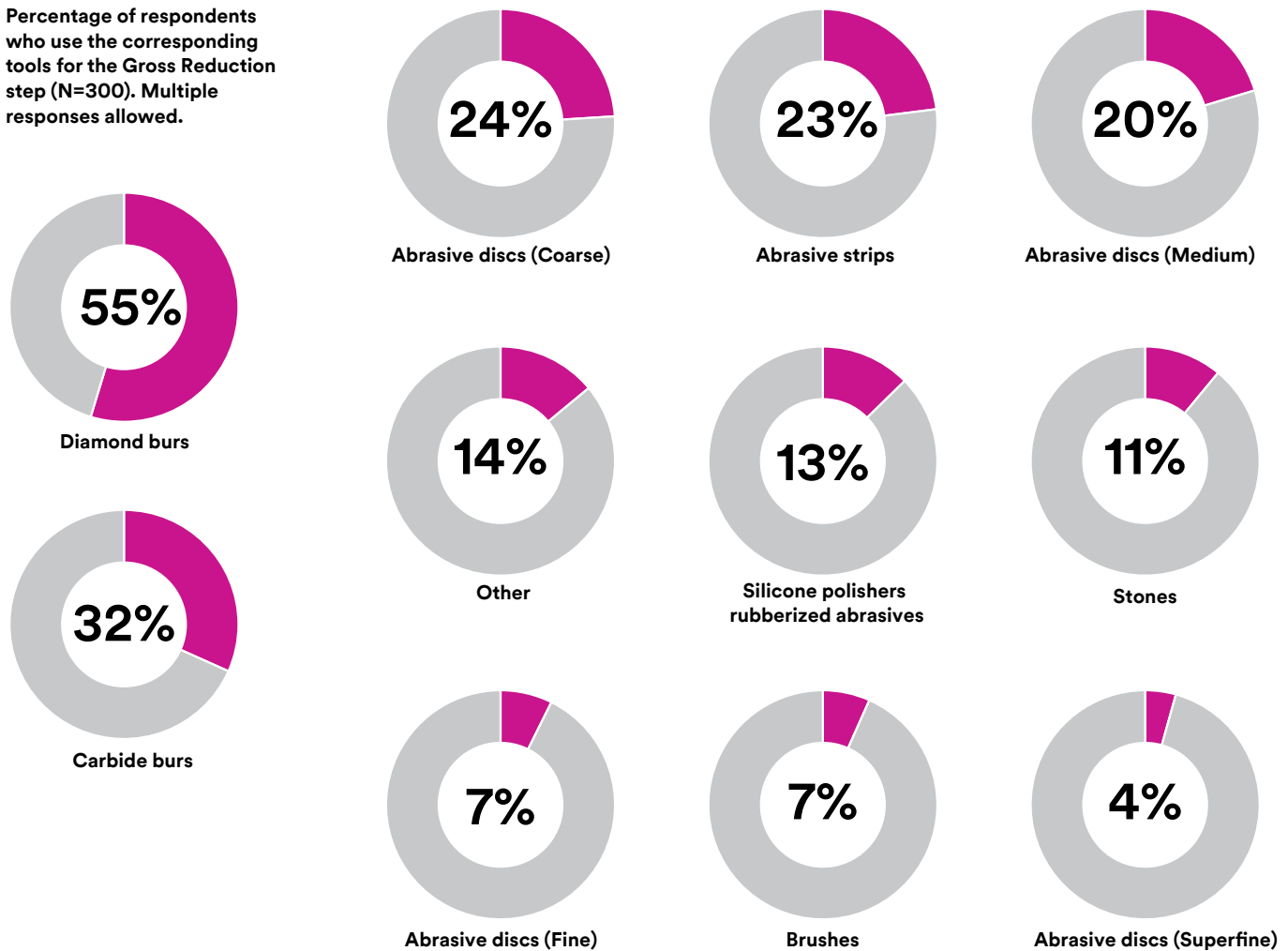


Contour:

The purpose of this step is to refine contours (size, shape, grooves, etc.) and margins (generating a smooth transition between tooth and restoration), reestablish contact with adjacent teeth to a normal and functional form and reduce surface roughness. At the end of this step, the restoration should have its desired form and a smooth, clean surface.

The most popular tools used for Final Contour are outlined below. Strips are also used for proximal areas. These tools are not as aggressive as those used for Gross Reduction, e.g., medium grit, because only small amounts of material are being removed from the restoration.

Percentage of respondents who use the corresponding tools for the Final Contour step (N=300). Multiple responses allowed.

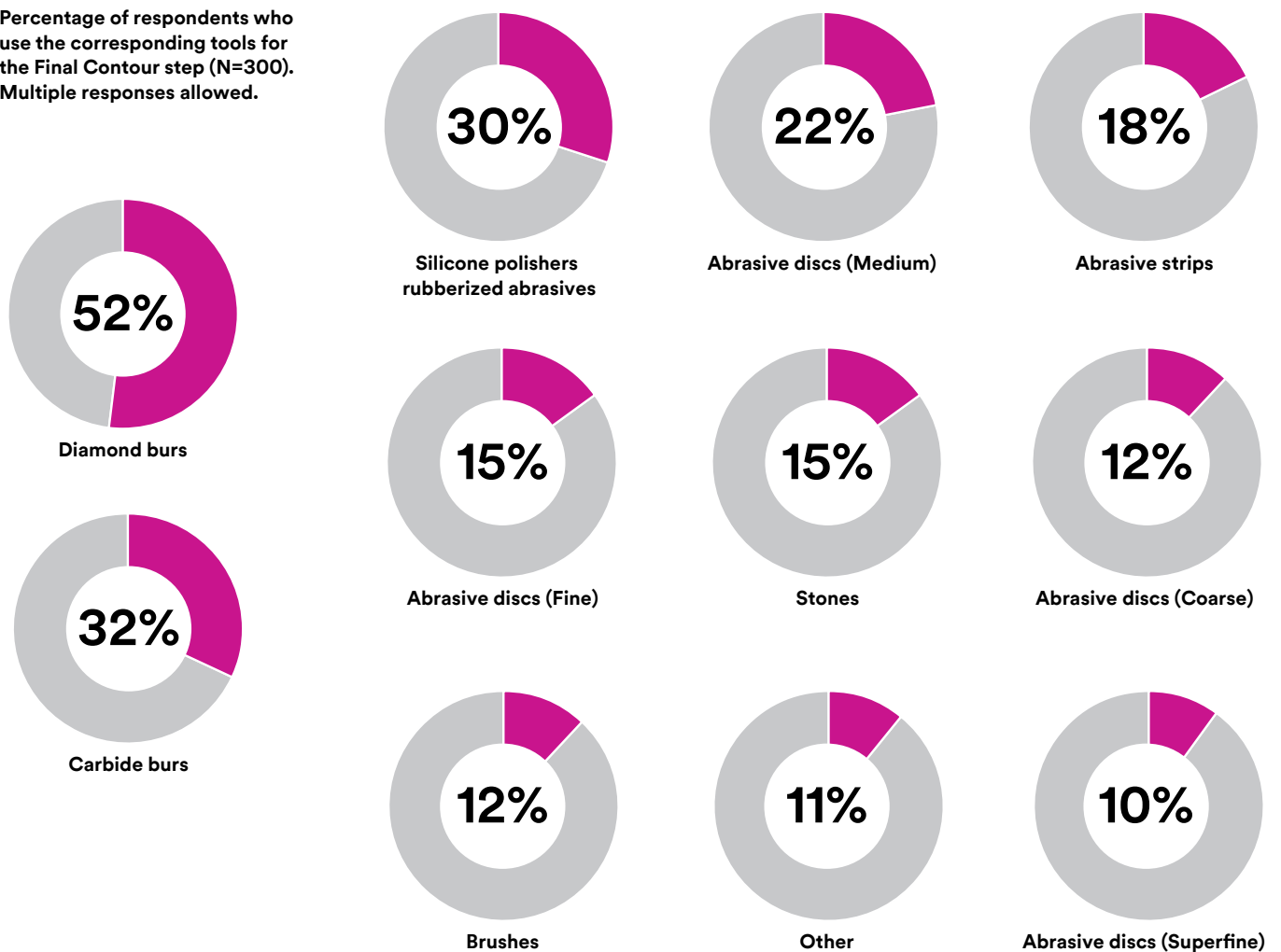


Finish: Smooth Surface (scratch removal)

This step reduces the scratch depth and/or removes lighter scratches produced during the Gross Reduction and Final Contour steps by the more aggressive tools. This step should leave the surface smoother and sometimes is considered a pre-polish step. The tools used for finishing (see below) are less aggressive than the previous step (e.g., a fine grit). Abrasive strips are also used for proximal areas. The Sof-Lex™ Pre-Polishing Spiral is designed to remove scratches and smooth the surface of teeth.

Responses from research—which instrument do you use to smooth?

Percentage of respondents who use the corresponding tools for the Pre-Polishing step (N=300). Multiple responses allowed.

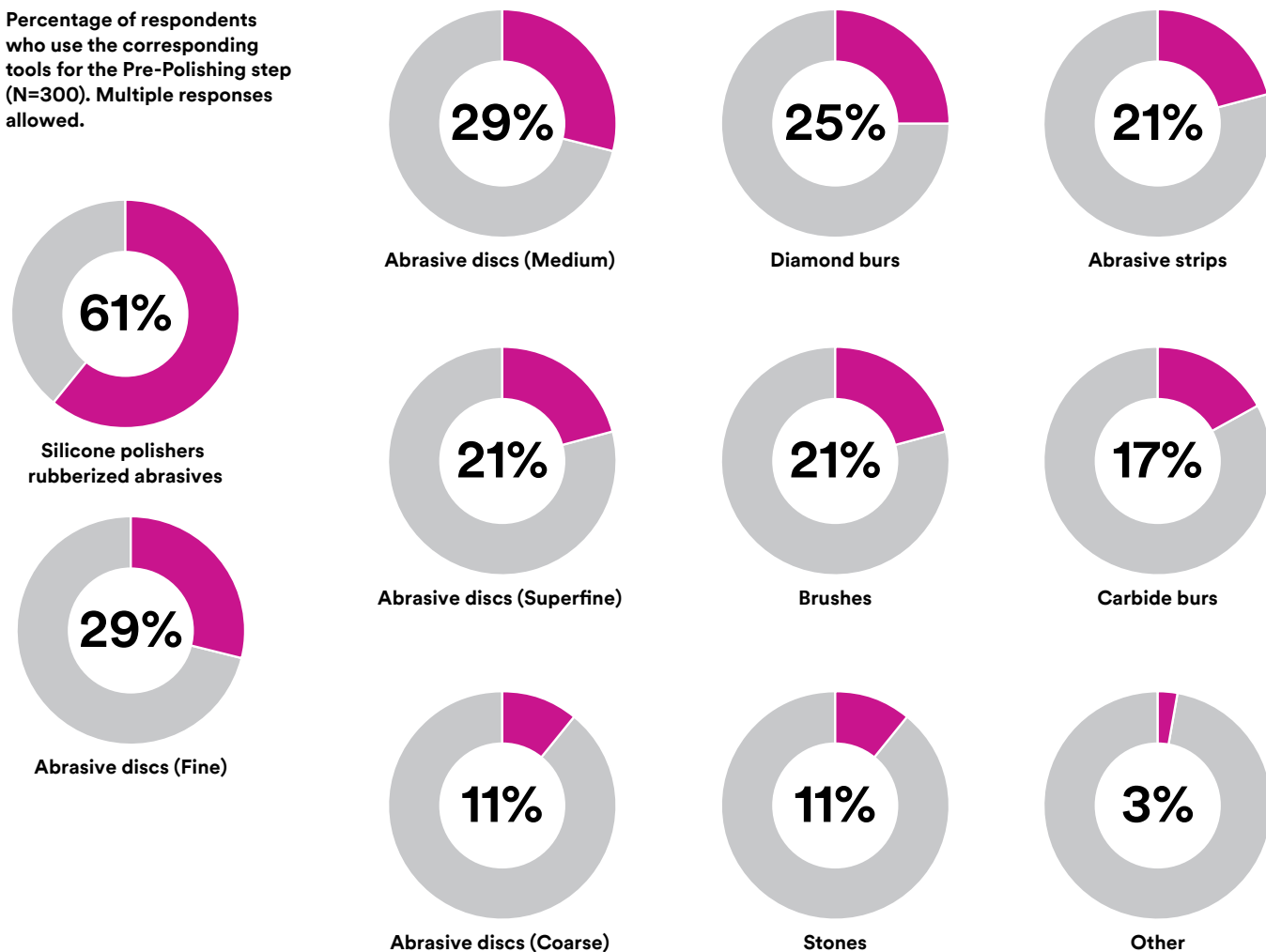

High-Gloss Polish

The objective of this step is to further smooth the surface to produce a high-gloss shine on the restoration. The tools used for this step are the least aggressive, e.g., super or ultra-fine grits. Strips are also used for proximal areas. The Sof-Lex Diamond Polishing Spiral produces a high-gloss polish.

Responses from research—which instrument do you use to polish?

Percentage of respondents who use the corresponding tools for the final High-Gloss Polish step (N=300). Multiple responses allowed.




Technology

The unique “bristle” shape is adapted for 3M Oral Care from the patented radial bristle disc design developed by the 3M Abrasives Division. This technology allows for:

- a continuous supply of mineral to work surfaces
- a variety of grits, diameters and thicknesses
- a flexible instrument that conforms to varied surfaces
- the instrument to generate minimal heat during use

3M Oral Care optimized this design for dental applications. The two abrasive spirals:

- Adapt to all tooth surfaces, eliminating the use of multiple tools (shapes) for a single purpose—to fit various contours
- Can be used on anterior and posterior restorations
- Beige—smoothes and removes scratches
- Pink—polishes the restoration to a high diamond paste-like gloss for a natural-looking surface
- Utilize the easy-to-use “Pop-On” 3M™ ESPE™ Mandrel
- Can be sterilized and reused multiple times

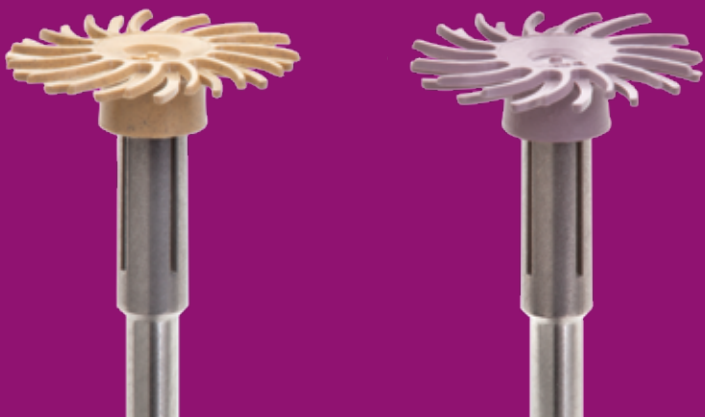

Performance

Performance for finishing and polishing systems are typically evaluated by measuring surface roughness, gloss measurement and surface appearance.

Surface Roughness

Surface roughness measures the smoothness of a surface. Typically, as a clinician uses a finisher and polisher sequence, the embedded abrasive particles become smaller in size with each step. The finer the abrasive particle, the smoother the surface should appear.

Gloss Measurements

Gloss is the shine or luster on a smooth surface. A smooth, uniform composite surface will reflect the greatest amount of light. A high gloss measurement indicates a shiny, more reflective surface.

Toothbrush abrasion is used as a method to simulate wear on any given composite. This test is used to measure how well a composite will retain polish.

Surface Appearance

Surface appearance is also observed to visualize how the abrasive particles are interacting with a composite surface. Gross reduction and contouring shape a composite restoration, but can cause deeper scratches to appear in the restoration surface. After smoothing and polishing, the surface should appear smooth and uniform in texture. Lack of a smooth surface can cause a restoration to appear dull or lacking in shine or gloss.

Sample Preparation


Samples of 3M™ ESPE™ Filtek™ Supreme Ultra Universal Composite and other brands of composite were prepared and then polished with 320 grit sandpaper. This uniform surface simulated a typical clinically relevant finished surface prior to polishing. The composite samples were then polished with various systems according to the manufacturer's instructions. Samples were polished with the 3M™ ESPE™ Sof-Lex™ Diamond Polishing System or the same brand polishing system as the composite tested.

Surface Roughness Measurement

Surface roughness can be quantified using surface profilometers. A profilometer drags a stylus over a composite surface and Ra (average surface roughness) values are recorded. Ra is the average surface roughness expressed in units of height over a set distance. The lower the Ra value, the smoother the surface.

In the chart below, samples of Filtek™ Supreme Ultra Universal Restorative were polished with a variety of polishing systems. Each instrument in the polishing system was used for 15 seconds and used according to the manufacturer’s instructions.

Surface Roughness (Ra) on Filtek™ Supreme Ultra Universal Restorative: Sof-Lex™ Diamond Polishing System vs. Other Polishing Systems

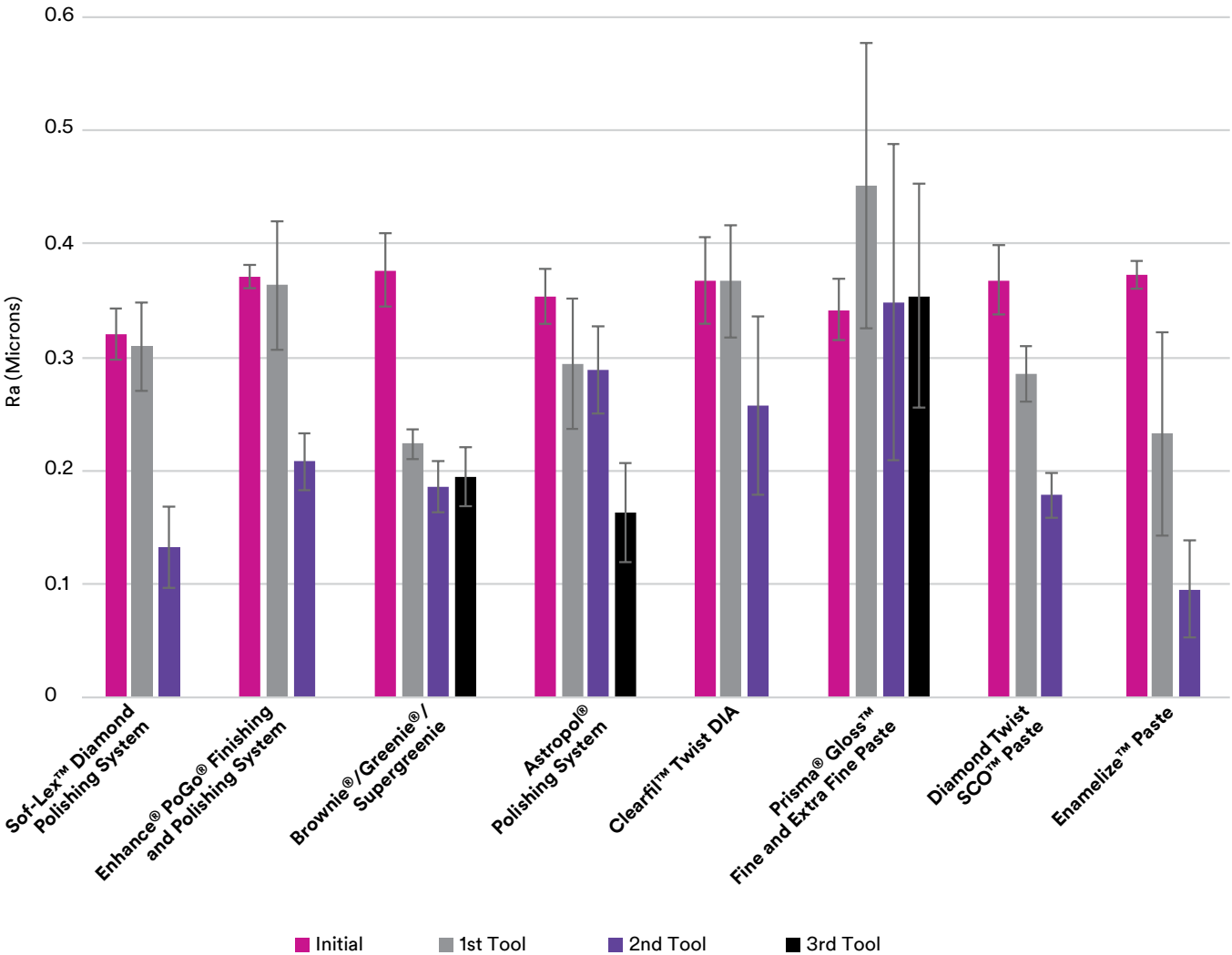


Gloss Measurements


Surface gloss or reflectance indicates how polished or shiny the surface can become after treatment. Gloss is measured using a gloss meter, which projects a beam of light at a fixed intensity and angle and measures the amount of reflected light at an equal but opposite angle.

In the first chart below, samples were prepared as previously described using Filtek™ Supreme Ultra Universal Restorative and polished with a variety of finishing and polishing systems. In the second chart, samples were prepared as previously described using a variety of composites and polished with Sof-Lex™ Diamond Polishing System. After final gloss measurements, toothbrush abrasion testing was also conducted to measure the polish retention of the particular composite.

Final Gloss on Filtek™ Supreme Ultra Universal Restorative: Sof-Lex™ Diamond Polishing System vs. Other Polishing Systems

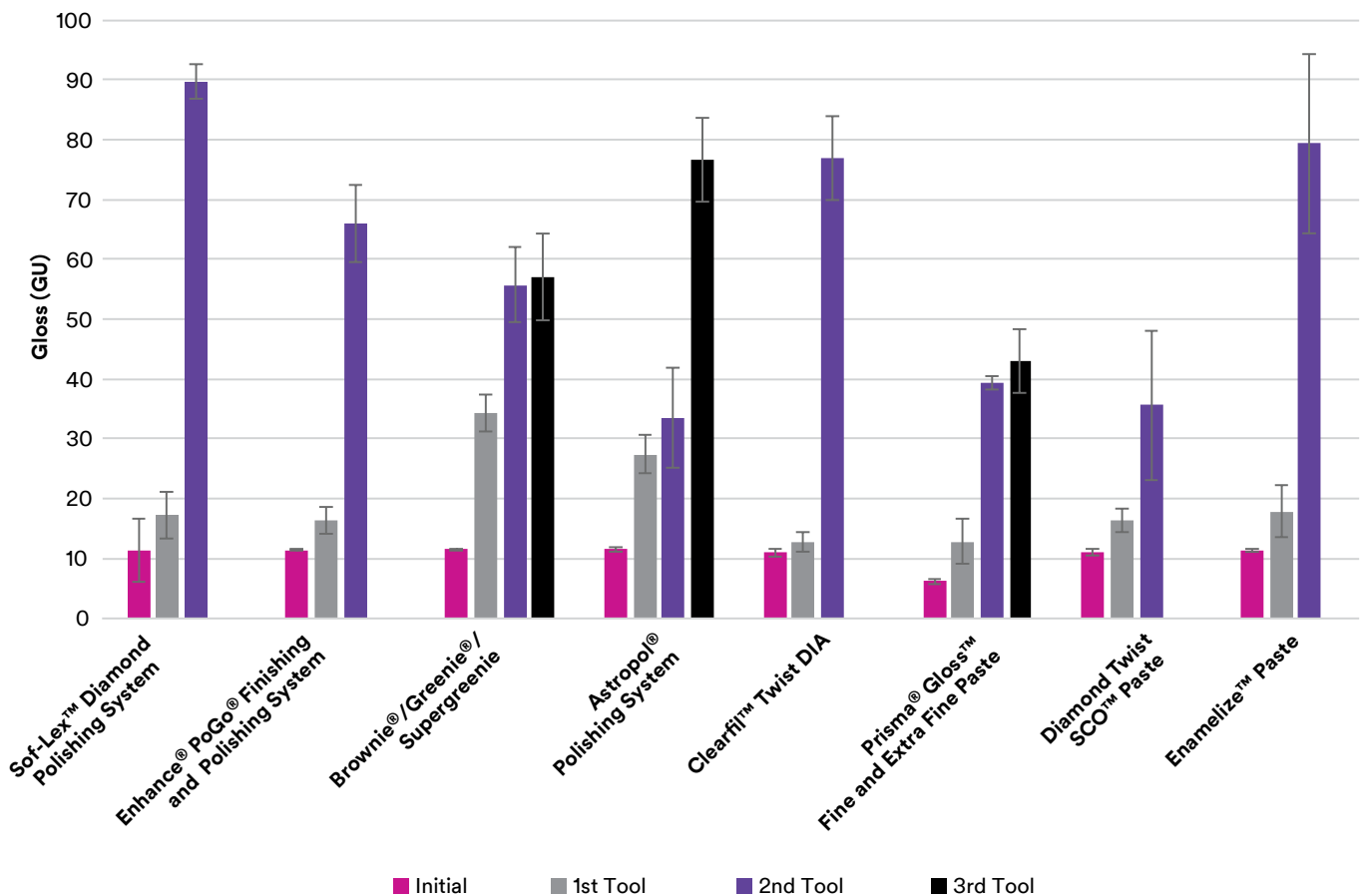

Final Gloss Before & After Toothbrush Abrasion: Sof-Lex™ Diamond Polishing System on Other Composites

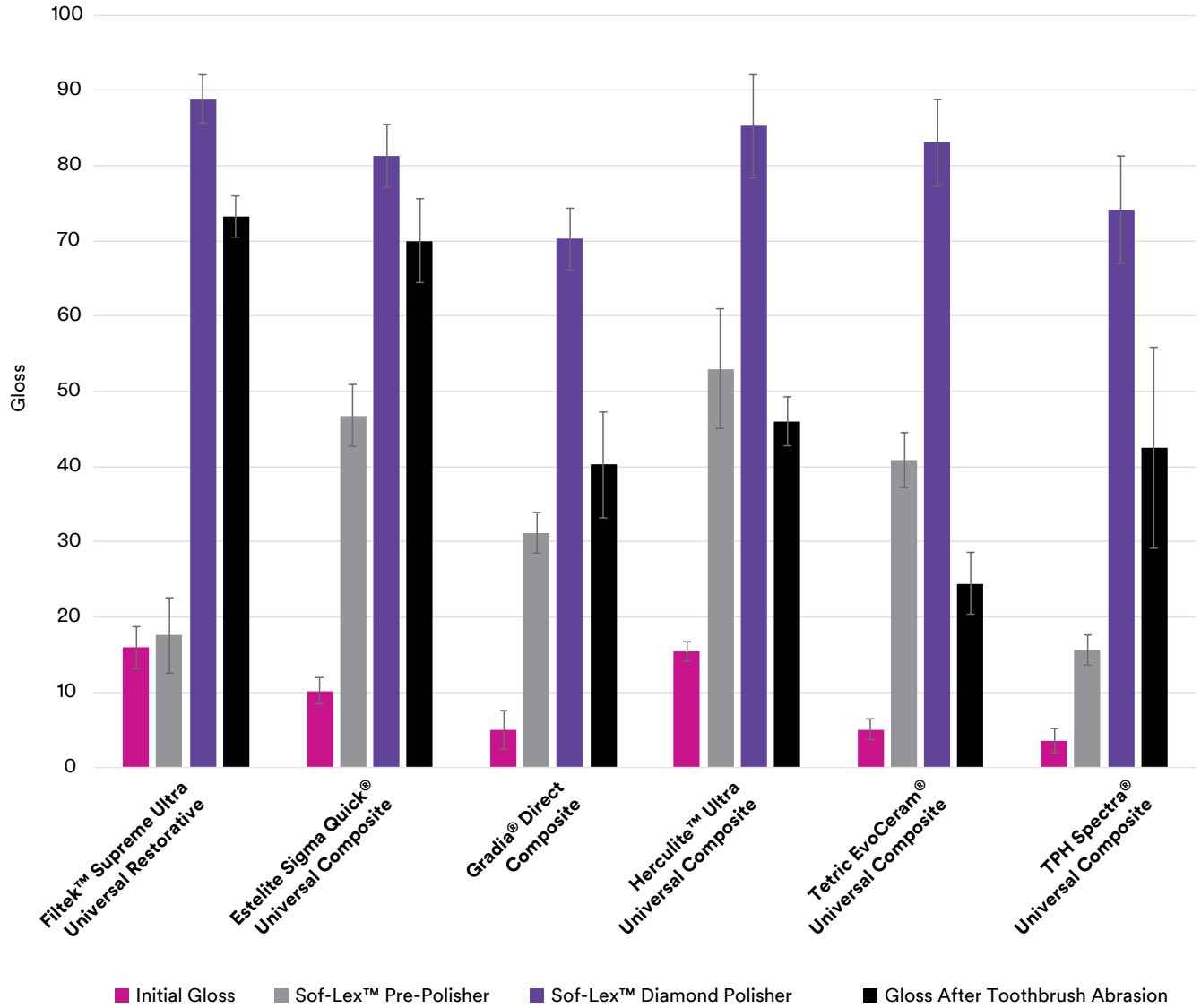

System Final Gloss

The following tiles visually show the difference in final gloss for each system. A more defined reflection suggests a higher final gloss.



Sof-Lex™ Diamond Polishing System on Filtek™ Supreme Ultra Universal Restorative (2 tools, 15 sec. each)



Occlubrush® Polishing System on Herculite™ Ultra Universal Composite (1 tool, 2 15 sec. steps)



Enhance® PoGo® Finishing and Polishing System on TPH Spectra® Universal Composite (2 tools, 15 sec. each)




Astropol® Polishing System on Tetric EvoCeram® Universal Composite (3 tools, 15 sec. each)

System Gloss Retention

Composite restorations may lose their gloss and show signs of wear over time.

The following composites were prepared using the same protocol as before. Samples were then polished with the corresponding polishing system (same manufacturer). After final polishing, samples received 6,000 cycles of toothbrush abrasion. Final gloss measurements were recorded.

Polish Retention: Samples Polished with the Corresponding Polishing System

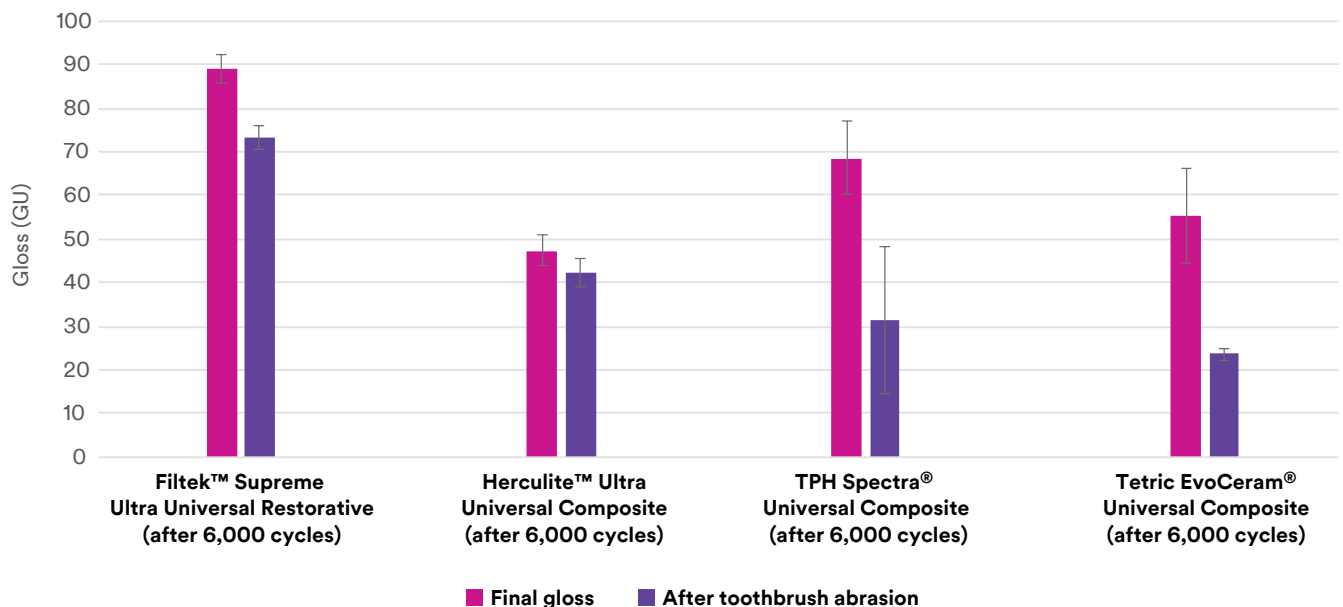


System Surface Roughness

In the chart below, a variety of composite samples were polished with corresponding polishing systems. Each instrument in the polishing system was used for 15 seconds according to the manufacturer’s instructions. The graph shows final surface roughness for each system.

Final Surface Roughness:

Composite with Same Brand Finishing & Polishing System

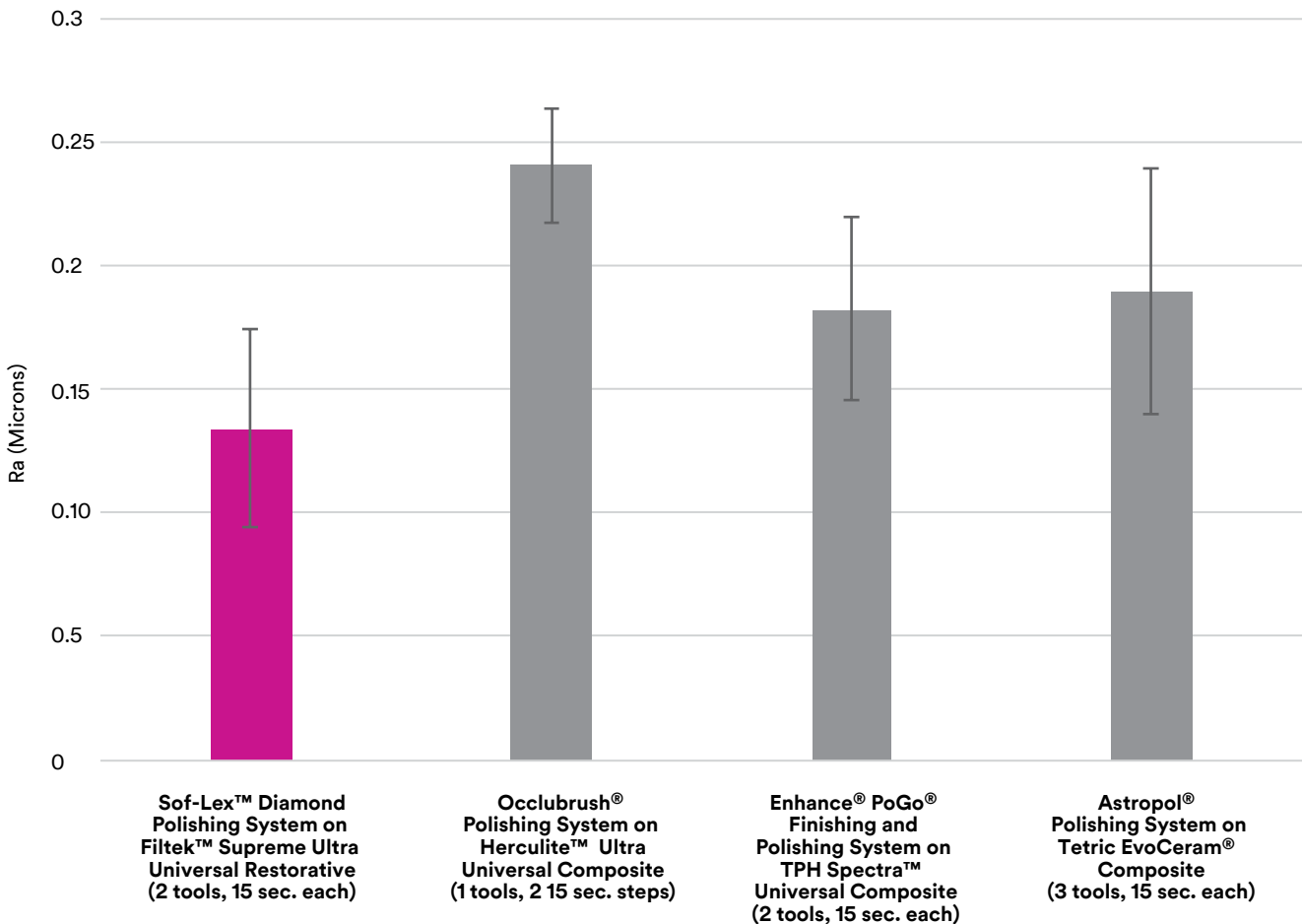


System Surface Visualization


Visual inspection of a composite can reveal how uniform and smooth a surface is after polishing. These polished surfaces are visualized using a scanning electron microscope (SEM). The TOPO setting is ideal for visualizing surface tomography is at a magnification of 400X.

Samples of Filtek™ Supreme Ultra Universal Restorative and a variety of other composites were prepared.


These samples were then polished with corresponding polishing systems following the manufacturer's instructions. SEM photos were captured to help visualize the surface.

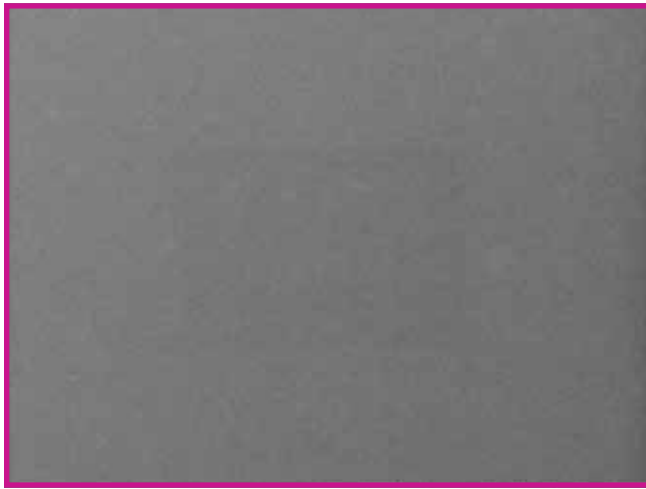


Filtek™ Supreme Ultra Universal Restorative polished with Sof-Lex™ Diamond Polishing System

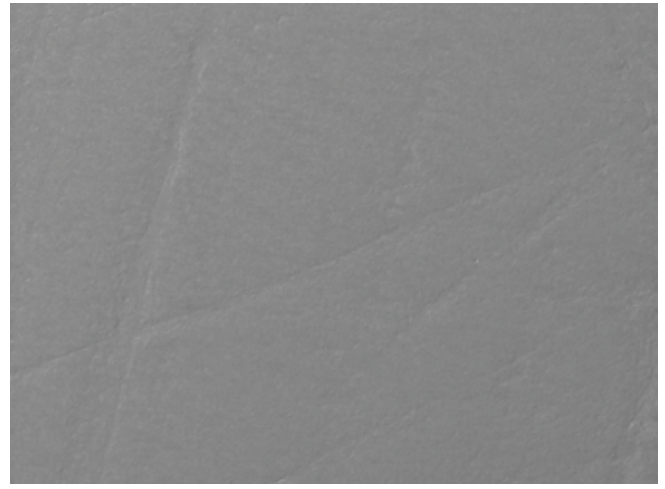

TPH Spectra® Universal Composite polished with Enhance® PoGo® Finishing and Polishing System

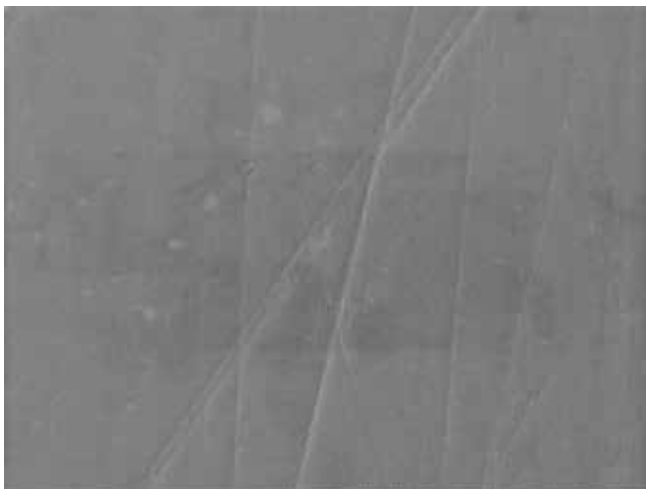

Tetric EvoCeram® Universal Composite polished with Astropol® Polishing System

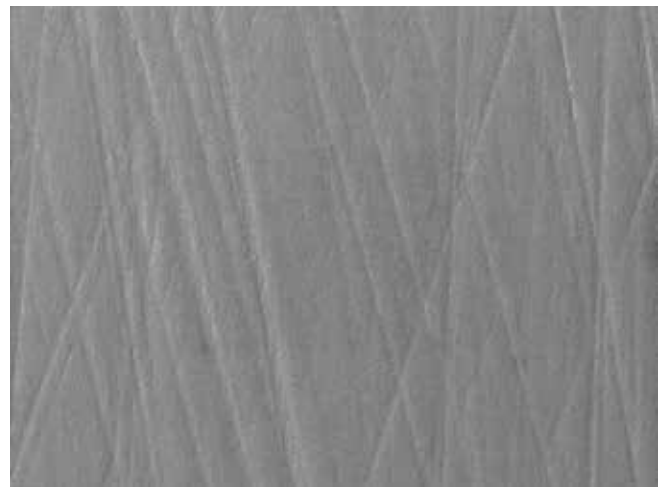

Herculite™ Ultra Universal Composite polished with Occlubrush® Polishing System

Multi-use

The Sof-Lex Diamond Polishing system can be sterilized and reused multiple times.

The well-known brands of rubberized polishers are single-use and must be disposed of after a single patient.

Cleaning and sterilization were validated and instructions are found in the instructions for use.


Ordering Information



Item #	Product Information
5092-I	Sof-Lex™ Diamond Polishing System—Introductory Kit Kit includes: 5 Sof-Lex™ Pre-Polishing Spirals (Beige); 5 Sof-Lex™ Diamond Polishing Spirals (Pink); Instructions for Use; Technique Card
5090	Sof-Lex™ Pre-Polishing Spirals—Refill Includes: 15 Sof-Lex™ Pre-Polishing Spirals (Beige)
5091	Sof-Lex™ Diamond Polishing Spirals—Refill Includes: 15 Sof-Lex™ Diamond Polishing Spirals (Pink)

Sof-Lex™ Contouring and Polishing Discs Refills


Each refill contains 85 discs.

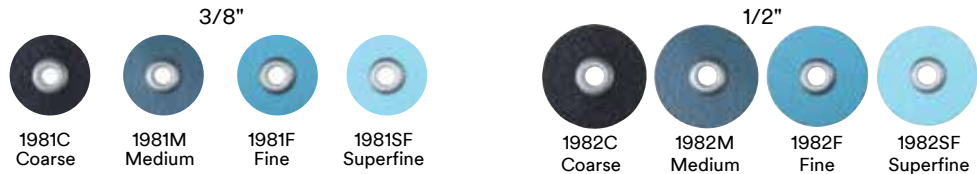


2380	Sof-Lex™ Extra-Thin Discs & Strips—System Kit Kit includes: 30 discs each of coarse, medium, fine and superfine grits in both 1/2" and 3/8" diameters; 1 RA mandrel.
------	---

Sof-Lex™ Extra-Thin Contouring and Polishing Discs Refills


Each refill contains 85 discs.

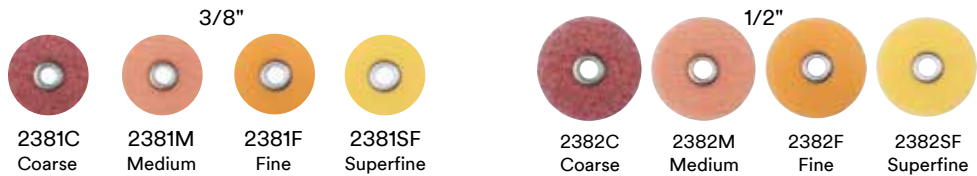


Sof-Lex™ Finishing Strips Refills	
1954	7" x 5/32" Clear Coarse/Medium Contains 150 strips

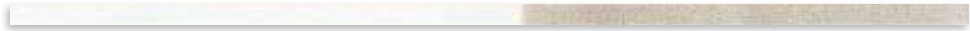

1954N	7" x 5/64" Clear Coarse/Medium—Narrow Contains 100 strips
-------	---

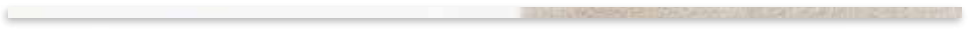

Sof-Lex™ Polishing Strip Refills	
1956	7" x 5/32" Yellow Fine/Superfine Contains 120 strips

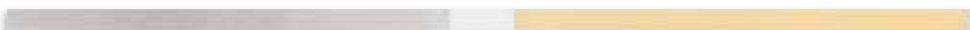

Mandrel Refills	
1983FG	3M™ ESPE™ FG Mandrels Refill (3)
1983RA	3M™ ESPE™ RA Mandrels Refill (3)
1983HP	3M™ ESPE™ HP Mandrels Refill (3)

3M.com/SofLex



3M Oral Care
2510 Conway Avenue
St. Paul, MN 55144-1000 USA
Phone 1-800-634-2249
Web 3M.com/dental

3M Canada
Post Office Box 5757
London, Ontario N6A 4T1
Canada
Phone 1-888-363-3685

3M, ESPE, Elipar and Sof-Lex are trademarks of 3M or 3M Deutschland GmbH. Used under license in Canada. Printed in USA. © 3M 2016. All rights reserved. 70-2013-0723-1